# Pocket BVM User Manual

## Adult Size, 1000 mL

### TABLE OF CONTENTS

|                                 |     |    |
|---------------------------------|-----|----|
| 1. English.....                 | 1   | EN |
| 2. French – Français.....       | 10  | FR |
| 3. German – Deutsch.....        | 20  | DE |
| 4. Spanish – Español.....       | 31  | ES |
| 5. Dutch – Nederlands.....      | 41  | NL |
| 6. Italian – Italiano.....      | 53  | IT |
| 7. Polish – Polski.....         | 63  | PL |
| 8. Danish – Dansk.....          | 74  | DA |
| 9. Norwegian – Norsk.....       | 84  | NO |
| 10. Czech – Čeština.....        | 94  | CS |
| 11. Swedish – Svenska.....      | 104 | SV |
| 12. Finnish – Suomi.....        | 114 | FI |
| 13. Icelandic – Íslenska.....   | 124 | IS |
| 14. Portuguese – Português..... | 134 | PT |

 **MANUFACTURED BY:**  
Micro BVM™ Systems Ltd.  
11 Ruzin St. Jerusalem 93870 Israel

 +972-2-6517191

 [info@microbvm.com](mailto:info@microbvm.com)

 [www.microbvm.com](http://www.microbvm.com)

MD

|                                                                                   |                                                                                                               |                |                      |
|-----------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------|----------------|----------------------|
|  | <b>Instruction for Use (IFU) – Pocket BVM™ manual resuscitator system</b><br><b>1000 ml: PBVM-KN, PBVM-KG</b> |                |                      |
|                                                                                   | <b>Doc # IFU-1000-EN</b>                                                                                      | <b>Rev. 01</b> | <b>Page 1 of 143</b> |

## 1. DEVICE DESCRIPTION

### Adult single patient use resuscitator including face mask, reservoir bag and oxygen inlet.

The Pocket BVM™ is a portable, operator-powered manual resuscitator. It is supplied as a complete set including a silicone face mask, a self-inflating resuscitator bag, an oxygen reservoir bag, an oxygen inlet connection and a 2 meter extension tube for oxygen supply tubing. The device is designed to be folded into a compact, robust polycarbonate case for easy storage and rapid deployment in emergency situations.

- **Sterility Status: The Pocket BVM is Non-Sterile.**
- **Service & Reuse:** The Pocket BVM is a single patient – single use device and, as such, no service is required. It should not be re-used on a same patient or others.

The device is compatible with the accessory Pocket BVM™ Extension tube of 75 cm length, which is not provided with the device.

## 2. INDICATIONS FOR USE & INTENDED USE

The Pocket BVM™ manual resuscitator is intended for manual pulmonary resuscitation and emergency respiratory support of patients by manually applying positive pressure to the airway to support ventilation and oxygenation. It acts as a temporary means of providing essential gas exchange until more definitive medical treatment can be initiated.

## 3. INTENDED USER PROFILE

The device is for use by medical personnel, emergency responders, and individuals specifically trained in manual ventilation and resuscitation techniques.

## 4. TARGET PATIENT POPULATION

This device is designed for adult patients with an ideal body mass **greater than 40 kg (88 lbs)**.

## 5. CONTRAINDICATIONS

The device is **contraindicated** for use on infants and children weighing less than 40 kg.

## 6. CLINICAL BENEFITS

The Pocket BVM provides emergency respiratory support, enabling effective artificial ventilation and oxygenation for patients experiencing respiratory arrest or insufficiency, thereby supporting vital gas exchange and reducing the risk of hypoxia.

|                                                                                   |                                                                                                               |                |                      |
|-----------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------|----------------|----------------------|
|  | <b>Instruction for Use (IFU) – Pocket BVM™ manual resuscitator system</b><br><b>1000 ml: PBVM-KN, PBVM-KG</b> |                |                      |
|                                                                                   | <b>Doc # IFU-1000-EN</b>                                                                                      | <b>Rev. 01</b> | <b>Page 2 of 143</b> |

## 7. WARNINGS

- **WARNING:** Use the correct size resuscitator for the ideal body mass of the patient (>40 kg) to avoid the risk of hypoventilation or barotrauma.
- **WARNING:** To avoid the risk of barotrauma, do not override the mechanical pressure relief unless it is clinically justified. Care needs to be taken to immediately restore the pressure relief function after the clinical need is resolved.
- **WARNING:** To avoid the risk of barotrauma, do not use abrupt and forceful compressions unless it is clinically justified because they can cause high airway pressures.
- **WARNING: This manual resuscitator is not equipped with a pressure-limiting valve.** The operator must closely monitor airway pressure and observe the patient's chest movement to prevent excessive pressure delivery, which could result in pulmonary barotrauma.
- **WARNING:** Open flames during resuscitation with oxygen are dangerous and are likely to result in fire or death. Do not allow open flames or sparks within 2 meters of the resuscitator or any oxygen-carrying accessories.
- **WARNING:** Do not lubricate fittings, connections, tubing, or other accessories of the resuscitator to avoid the risk of fire and burns.
- **WARNING:** Reduced or insufficient airflow may result in brain damage to the patient being ventilated.
- **WARNING:** Improper function could lead to inadequate ventilation, resulting in serious injury or death.
- **WARNING:** If efficient ventilation cannot be obtained with this device, the operator should immediately replace to alternative device or switch to mouth-to-mouth and/or mouth-to-nose resuscitation.
- **WARNING:** Use in contaminated environments can be hazardous unless entrainment of contaminants is prevented.
- **WARNING:** Patient expired gas is potentially infectious. Breathing filters can reduce but not eliminate contamination risk.
- **WARNING:** Incorrect use of this device can lead to excessive patient rebreathing and death.
- **WARNING:** Oil or grease should not be used in close proximity to oxygen equipment. Fire may result.

|                                                                                   |                                                                                                               |                |                      |
|-----------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------|----------------|----------------------|
|  | <b>Instruction for Use (IFU) – Pocket BVM™ manual resuscitator system</b><br><b>1000 ml: PBVM-KN, PBVM-KG</b> |                |                      |
|                                                                                   | <b>Doc # IFU-1000-EN</b>                                                                                      | <b>Rev. 01</b> | <b>Page 3 of 143</b> |

- **WARNING:** Avoid using an oxygen concentration more than that which is clinically required by the patient. Delivering excessive oxygen can increase the risk of oxygen toxicity.
- **WARNING:** Check all accessory connections before use for a secure fit before use to avoid potential disconnection during operation.
- **WARNING:** Do not use with other devices any of the parts or accessories provided with the Pocket BVM™ device.
- **WARNING:** Do not pull or apply excessive force to the resuscitator bag or the patient valve housing during use, as this may compromise connection stability with accessories.


## 8. PRECAUTIONS & CAUTIONS

- **CAUTION:** The Pocket BVM is intended for Single Patient Use only. It must not be re-processed, cleaned, or reused. Reuse of the device may lead to cross-contamination and mechanical failure, compromising patient safety.
- **CAUTION:** Never store the resuscitator in a squeezed or compressed state, other than folded at the folding zones in the original packaging. Permanent deformation of the bag will occur, which may reduce the ventilation efficiency. Always store the resuscitator in its original box.
- **CAUTION:** Avoid contamination with oil and grease. Such products may affect the integrity of the resuscitator materials.
- **CAUTION:** Do not use the device if the original protective packaging is damaged or if the expiration date has passed.
- **CAUTION:** U.S. Federal law restricts this device to sale by or on the order of a licensed healthcare practitioner.

|                                                                                   |                                                                                                               |                |                      |
|-----------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------|----------------|----------------------|
|  | <b>Instruction for Use (IFU) – Pocket BVM™ manual resuscitator system</b><br><b>1000 ml: PBVM-KN, PBVM-KG</b> |                |                      |
|                                                                                   | <b>Doc # IFU-1000-EN</b>                                                                                      | <b>Rev. 01</b> | <b>Page 4 of 143</b> |

## 9. INSTRUCTIONS FOR USE

### 9.1 Unpack and Assemble the Pocket BVM (Figure 1)

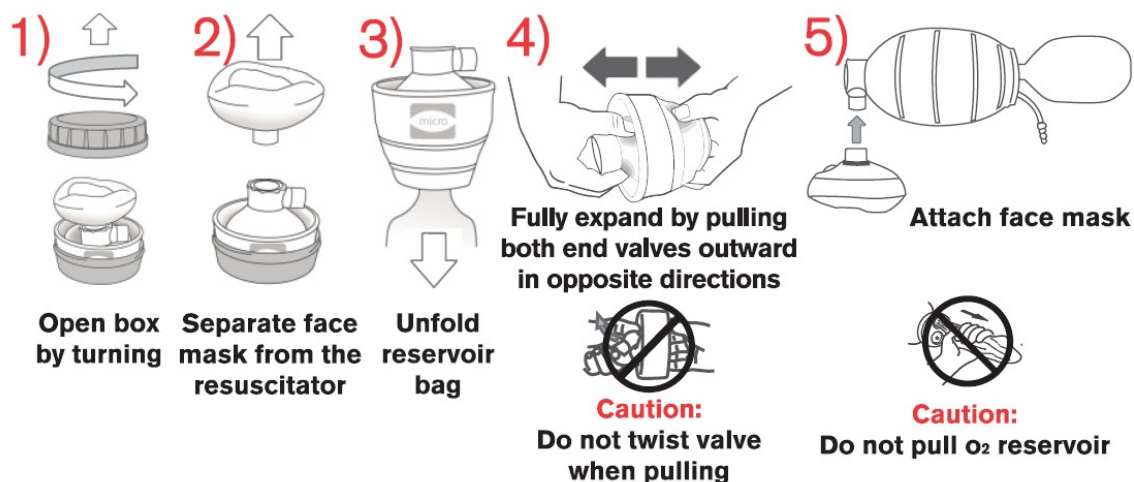**Figure 1**

1. Open the box by turning the cover.
  2. Separate the Face Mask from the folded Resuscitator Bag.
  3. Unfold the O<sub>2</sub> reservoir bag.
  4. Fully expand the Resuscitator Bag from its folded position by pulling both end valves outwards in opposite directions.
- Note:** Do not twist valves while pulling; do not pull the resuscitator bag or the patient valve housing to avoid accidental disconnection of parts.
5. Attach the Face Mask to the Patient Valve housing.
  6. Ensure a secure mechanical connection for all attached accessories prior to use; the integrated sealing radius is designed to provide a robust, airtight seal during manual operation.

### 9.2 Functional Test for the Pocket BVM, Prior to Use:

- Check the resuscitator by squeezing a few times to confirm that the Inlet Valve and the Patient Valve are both functioning (i.e., air is entering the bag and exiting from the patient valve).
- Check that the bag fills through visible expansion and relaxation.

|                                                                                   |                                                                                                               |                |                      |
|-----------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------|----------------|----------------------|
|  | <b>Instruction for Use (IFU) – Pocket BVM™ manual resuscitator system</b><br><b>1000 ml: PBVM-KN, PBVM-KG</b> |                |                      |
|                                                                                   | <b>Doc # IFU-1000-EN</b>                                                                                      | <b>Rev. 01</b> | <b>Page 5 of 143</b> |

- Block the Patient Port with your thumb and squeeze the bag and hold for positive pressure to verify the system is leak-free.
- If there is a leak from the 'mask filter housing', try to tighten the 'peep holder' part clockwise.
- If any components (such as the valves) are not fully functional, **do not use** this device.
- All accessories not supplied with the device must be verified for compatibility and ensure a secure connection prior to clinical use. Use of incompatible parts can result in degraded performance.

### 9.3 Ventilating the Patient:

- If resuscitating with Oxygen, attach the oxygen supply tube to a regulated oxygen source and set desired flow.

**Note:** *O<sub>2</sub> should be administered as soon as possible.*

- Use recommended techniques to position the patient correctly, to open and clear the mouth and airways, and to hold the mask firmly against the face.
- Ensure that the facemask forms an adequate seal with the patient's face. A tight seal is essential; failure may result in reduced or no ventilation.
- The Pocket BVM functions must be visible and audible: Listen for expiratory airflow from the valves to confirm function.
- During insufflation, observe the rise of the patient's chest. Release the bag abruptly and observe the lowering of the chest. If continued resistance is encountered, check for airway obstruction or correct backward head tilt.
- The correct ventilation frequency may vary: Please follow the current ventilation frequency recommended by regional, national or international guidelines.

### 9.4 Emergency Procedure (Vomit Contamination): In an emergency where the resuscitator is soiled with vomitus or secretions:

1. **Disconnect** the device from the patient.
2. **Clear** the mask and patient valve by shaking or wiping out the vomitus.
3. **Disassemble/Reassemble:** If necessary, briefly separate the patient valve from the bag to clear blockages, then promptly reassemble it.
4. **Functional Test:** Squeeze the bag once to confirm proper operation before reconnecting.

|                                                                                   |                                                                                                               |                |                      |
|-----------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------|----------------|----------------------|
|  | <b>Instruction for Use (IFU) – Pocket BVM™ manual resuscitator system</b><br><b>1000 ml: PBVM-KN, PBVM-KG</b> |                |                      |
|                                                                                   | <b>Doc # IFU-1000-EN</b>                                                                                      | <b>Rev. 01</b> | <b>Page 6 of 143</b> |

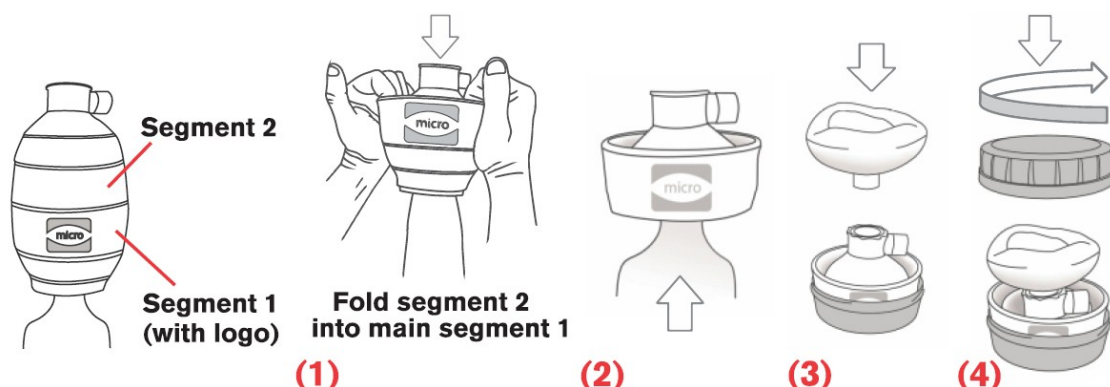
## 10. TECHNICAL SPECIFICATIONS

| Parameter                                                     | Specification Value                                                        |
|---------------------------------------------------------------|----------------------------------------------------------------------------|
| Body mass range for use                                       | >40 kg (88 lbs.)                                                           |
| Operating environmental limits                                | -18°C (0°F) to +50°C (+122°F); 15% to 85% RH                               |
| Storage environmental limits                                  | -40°C (-40°F) to +60°C (+140°F)                                            |
| Substances which can be administered                          | Oxygen                                                                     |
| Delivered Oxygen Concentration                                | 94.7% at an input flowrate ≤ 15l/min                                       |
| Connector Dimensions (nominal)                                | 15 mm ID / 22 mm OD                                                        |
| Gas inlet tube connection                                     | 15 mm length x 6 mm O.D.                                                   |
| Tidal volume ( $V_T$ ) one hand                               | 600 ml                                                                     |
| Resuscitator Deadspace                                        | 4.0 ml                                                                     |
| Expiratory resistance (at 50 l/min)                           | 3.6 cmH <sub>2</sub> O                                                     |
| Inspiratory resistance (at 50 l/min)                          | 4.3 cmH <sub>2</sub> O                                                     |
| Maximum limited pressure                                      | ~15 cmH <sub>2</sub> O (structural limit)                                  |
| Gas source excessive flow (PEEP)                              | <2.1 cmH <sub>2</sub> O at 30 l/min                                        |
| External dimensions (Deployed)                                | 191 mm (L) x 108 mm (D)                                                    |
| External dimensions (In Case)                                 | 63 mm (H) x 134 mm (D)                                                     |
| Mass of the resuscitator, face mask and the resuscitator case | Resuscitator 142 g, Case 100 g, Mask 95 g.<br>Total < 470 g                |
| Shelf Life                                                    | 5 Years                                                                    |
| Operational Lifetime                                          | The device is designed for at least 24 hours of continuous operational use |

**Note regarding Oxygen Configuration:** To achieve the maximum delivered oxygen concentration, the resuscitator must be used in its full configuration (with the Oxygen Reservoir Bag attached and connected to a regulated medical oxygen source of 90%-100% O<sub>2</sub>). Using an oxygen concentrator will proportionally reduce the delivered oxygen concentration.

|                                                                                   |                                                                                                               |                |                      |
|-----------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------|----------------|----------------------|
|  | <b>Instruction for Use (IFU) – Pocket BVM™ manual resuscitator system</b><br><b>1000 ml: PBVM-KN, PBVM-KG</b> |                |                      |
|                                                                                   | <b>Doc # IFU-1000-EN</b>                                                                                      | <b>Rev. 01</b> | <b>Page 7 of 143</b> |

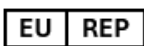
## 11. FOLDING AND PACKING FOR STORAGE

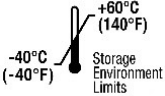
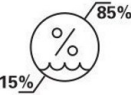
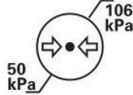
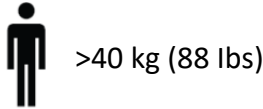
## 12. REPORTING SERIOUS INCIDENTS

Any serious incident that has occurred in relation to the device should be reported to the manufacturer at [info@microbvm.com](mailto:info@microbvm.com) and the competent authority of the Member State in which the user and/or patient is established.

## 13. SYMBOL GLOSSARY

| Symbol                                                                              | Title                                       | Description                                                                                |
|-------------------------------------------------------------------------------------|---------------------------------------------|--------------------------------------------------------------------------------------------|
|  | Manufacturer                                | Indicates the medical device manufacturer                                                  |
|  | Authorized Representative in European Union | Indicates the authorized representative in the European Union                              |
|  | Date of Manufacture                         | Indicates the date when the medical device was manufactured                                |
|  | Use-by date                                 | Indicates the date after which the medical device is not to be used                        |
|  | Batch Code                                  | Indicates the manufacturer's batch code so that the batch or lot can be identified         |
|  | Catalogue Number                            | Indicates the manufacturer's catalogue number so that the medical device can be identified |
|  | Importer                                    | Indicates the entity importing the medical device into the locale                          |
|  | Country of manufacture                      | To identify the country of manufacture of products                                         |

|                                                                                   |                                                                                                               |                |                      |
|-----------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------|----------------|----------------------|
|  | <b>Instruction for Use (IFU) – Pocket BVM™ manual resuscitator system</b><br><b>1000 ml: PBVM-KN, PBVM-KG</b> |                |                      |
|                                                                                   | <b>Doc # IFU-1000-EN</b>                                                                                      | <b>Rev. 01</b> | <b>Page 8 of 143</b> |

| Symbol                                                                              | Title                                                                  | Description                                                                                                                                                                       |
|-------------------------------------------------------------------------------------|------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|    | Non-sterile                                                            | Indicates a medical device that has not been subjected to a sterilization process                                                                                                 |
|    | Do not use if package is damaged and consult instructions for use      | Indicates a medical device that should not be used if the package has been damaged or opened and that the user should consult the instructions for use for additional information |
|    | Fragile, handle with care                                              | Indicates a medical device that can be broken or damaged if not handled carefully                                                                                                 |
|    | Keep away from sunlight                                                | Indicates a medical device that needs protection from light sources                                                                                                               |
|   | Keep dry                                                               | Indicates a medical device that needs to be protected from moisture                                                                                                               |
|  | Temperature Limit                                                      | Indicates the temperature limits to which the medical device can be safely exposed                                                                                                |
|  | Humidity limitation                                                    | Indicates the range of humidity to which the medical device can be safely exposed                                                                                                 |
|  | Atmospheric pressure limitation                                        | Indicates the range of atmospheric pressure to which the medical device can be safely exposed                                                                                     |
|  | Do not re-use                                                          | Indicates a medical device that is intended for one single use only                                                                                                               |
|  | Consult Instruction for use or consult electronic instructions for use | Indicates the need for the user to consult the instructions for use. System requirements to view the eIFU: Internet access and PDF reader software                                |
|  | Does not contain natural rubber latex                                  | Indicates that natural rubber latex was not used in the manufacturing of the product, its container, or its packaging                                                             |
|  | Body mass range for use                                                | Indicates the body mass range for use > 40 kg (88 lbs)                                                                                                                            |

|                                                                                   |                                                                                                               |                |                      |
|-----------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------|----------------|----------------------|
|  | <b>Instruction for Use (IFU) – Pocket BVM™ manual resuscitator system</b><br><b>1000 ml: PBVM-KN, PBVM-KG</b> |                |                      |
|                                                                                   | <b>Doc # IFU-1000-EN</b>                                                                                      | <b>Rev. 01</b> | <b>Page 9 of 143</b> |

| Symbol                                                                            | Title                          | Description                                                                                                      |
|-----------------------------------------------------------------------------------|--------------------------------|------------------------------------------------------------------------------------------------------------------|
|  | Medical Device                 | Indicates the item is a medical device                                                                           |
|  | Unique Device Identifier       | Indicates a carrier that contains Unique Device Identifier information                                           |
|  | Prescription Only              | Caution: Federal (USA) law restricts this device to sale by or on the order of a physician                       |
|  | Magnetic Resonance (MR) Unsafe | An item which poses unacceptable risks to the patient, medical staff or other persons within the MR environment. |
|  | CE Marking                     | Indicates that the medical device complies with the requirements of applicable European Union legislation        |

#### ACCESS TO ELECTRONIC INSTRUCTIONS (eIFU):

To view this document in other languages, scan the QR code below or visit

[www.microbvm.com](http://www.microbvm.com).

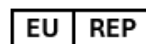
System requirements: Internet access and PDF reader software.

A paper copy of these instructions is available upon request at no cost and will be provided within 7 calendar days.

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